
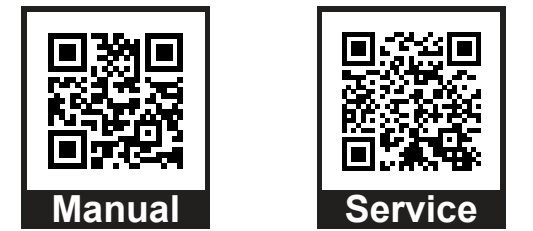
medisana®

MC 600

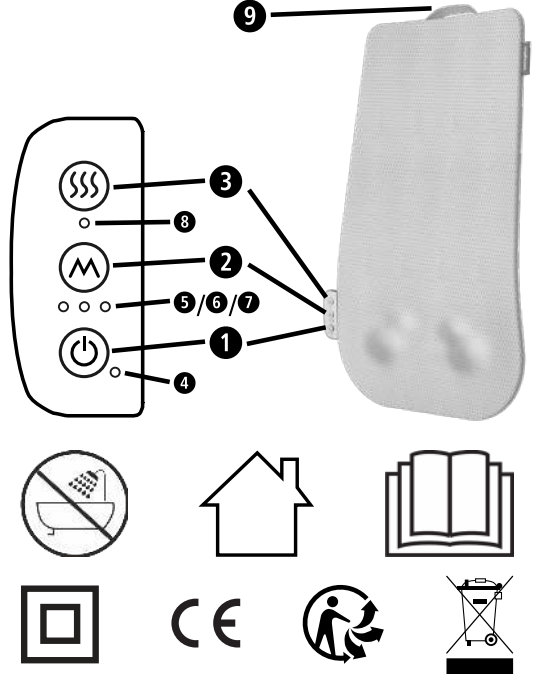medisana
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

REF 88967 **EAN** 40 15588 88967 7

Massagesitzauflage MC 600 **DE**

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Massage Seat Cover MC 600

- 1 On/Off button / 1 LED indicator for operation
- 2 Mode button (4 functions available)
3 3 LED indicator for function mode:
 - full back: 3 is lit
 - upper back: 3 is lit
 - lower back: 3 is lit
 - Spot massage: no LED is lit
- 4 Heat button / 2 LED indicator for heat
- 5 Carrying handle

Safety Instructions


Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- The device is not to be used for commercial use or medical purposes.
- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the power supply unit into a socket when the device is switched off.
- This device may only be operated with the mains adapter (XZ1200-2500G) supplied.
- Keep the power supply unit, cable and appliance well away from heat, hot surfaces, moisture and liquids.
- Do not touch the power supply unit if you are standing in water or if your hands are wet or damp.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Only plug the device into a socket where the plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power supply unit.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device by the cable or the power supply unit.
- You must not continue to use the unit if the cable or power supply unit are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that no one can trip over the cable. The cable must not be kinked, trapped or twisted.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Do not use this device to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the device, or only after consulting your doctor, if:
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Take care when using the heating function. The device has hot surfaces. Do not use the device on infants, vulnerable individuals or people who are not sensitive to heat.
- Before using the device, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.
- Check the power supply unit, cable and Massage Seat Cover carefully for damage before each use. Do not attempt to use a defective unit.
- Do not use the device if it or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not keep or place any electrical device on a surface where it is in danger of falling into a bathtub or sink.
- Only use the device in the proper place, as described in this instruction manual.
- Do not use the device while you are asleep or when you are lying in bed.
- Do not use the device before going to sleep. The massage has a stimulating effect. Never use the device while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Longer use shortens the life of the device and can lead to tension of the muscles instead of relaxation.
- Very long periods of use can lead to overheating. Always allow the device to cool down before using it again.
- Never store or use the device next to an electric stove or other sources of heat.
- Never leave the device unattended if it is connected to the power supply.
- Do not allow the device to come into contact with pointed or sharp objects.
- Do not cover the device when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
- Never sit or stand on the moving parts of the device with all your weight and do not place objects on the device.
- Do not use an external time switch or a separate remote control system to operate the device.
- Do not use the device whilst using aerosol sprays or simultaneously using oxygen.
- The device is maintenance free. Do not try to open it. There are no parts inside that are to be cleaned or serviced by the user.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all warranty rights. Consult your specialist dealer and have repairs carried out by the service centre to avoid dangers.
- Do not immerse the device in water or other fluids. If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
- Cleaning and user maintenance shall not be carried out by children without supervision.
- Before cleaning the unit, make sure it is switched off and disconnected from the mains. Clean the unit once it has cooled down using a soft, damp cloth.
- Under no circumstances should you use aggressive detergents, strong brushes, solvents or alcohol, as they may damage the surface of the unit.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, cool and dry place.

Intended use

The device is used to massage certain body parts for relaxation and well being. It is not intended for medical therapeutic use. It is not suitable for clinics, massage practices, wellness studios etc. The massage device is not a medical device and is not intended for the treatment of pain. In case of doubt or uncertainty, consult your doctor before use.

Voltage supply: via mains adapter, type XZ1200-2500G
Input: 100–240V~ 50/60Hz / Output: 12V $\pm$ 2.5A

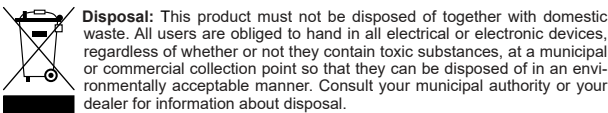Disposal: This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

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